


QUICKTRACH

Cricothyrotomy is an emergency procedure to avoid death from suffocation in case of upper airway obstruction if intubation or tracheotomy are impossible. VBM offers two Quicktrach-Sets with the following advantages:

- fast and safe airway in case of upper airway obstruction
- sterile sets are pre-assembled and immediately ready for use
- specially grinded needle tip makes prior incision via scalpel unnecessary and reduces the risk of bleeding
- stopper prevents needle from being inserted too deep to reduce the risk of posterior tracheal wall perforation
- necktape allows a fast and safe fixation
- quickest airway access via one-step procedure

Quicktrach II with Cuff

- thin cuff seals the trachea and allows efficient ventilation with aspiration protection
- stopper and safety clip reduce the risk of posterior tracheal wall perforation
- anatomically shaped cannula adjusts to the trachea due to “memory effect“

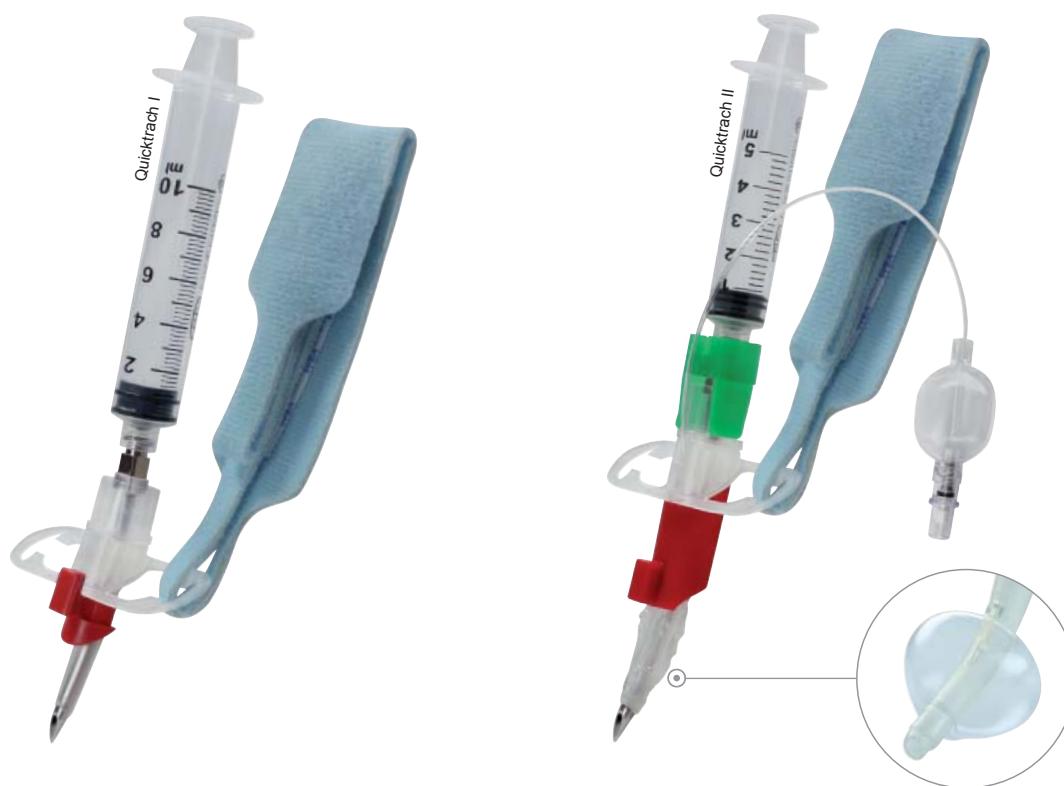

in blister packaging
for single use, sterile

Patient	I.D.	Quicktrach I	Quicktrach II with Cuff	Box
Infant	1.5 mm	REF 30-04-015-1	—	1
Child	2.0 mm	REF 30-04-002-1	—	1
Adult	4.0 mm	REF 30-04-004-1	REF 30-10-004-1	1



in plastic tube
for single use, sterile

- ideal for emergency bags (compact and robust)

Patient	I.D.	Quicktrach I	Quicktrach II with Cuff	Box
Child	2.0 mm	REF 30-04-902-1	—	1
Adult	4.0 mm	REF 30-04-904-1	REF 30-10-904-1	1

