


Cricothyrotomy

Cricothyrotomy-Set Quicktrach

Cricothyrotomy is an emergency procedure to avoid death from suffocation in case of upper airway obstruction if intubation or tracheotomy are impossible. VBM offers two Quicktrach-Sets with the following advantages:

- Fast and safe airway in case of upper airway obstruction
- Sterile sets are pre-assembled and immediately ready for use
- Specially grinded needle tip makes prior incision via scalpel unnecessary and reduces the risk of bleeding
- Stopper prevents needle from being inserted too deep to reduce the risk of posterior tracheal wall perforation
- Necktape allows a fast and safe fixation
- Quickest airway access via one-step procedure

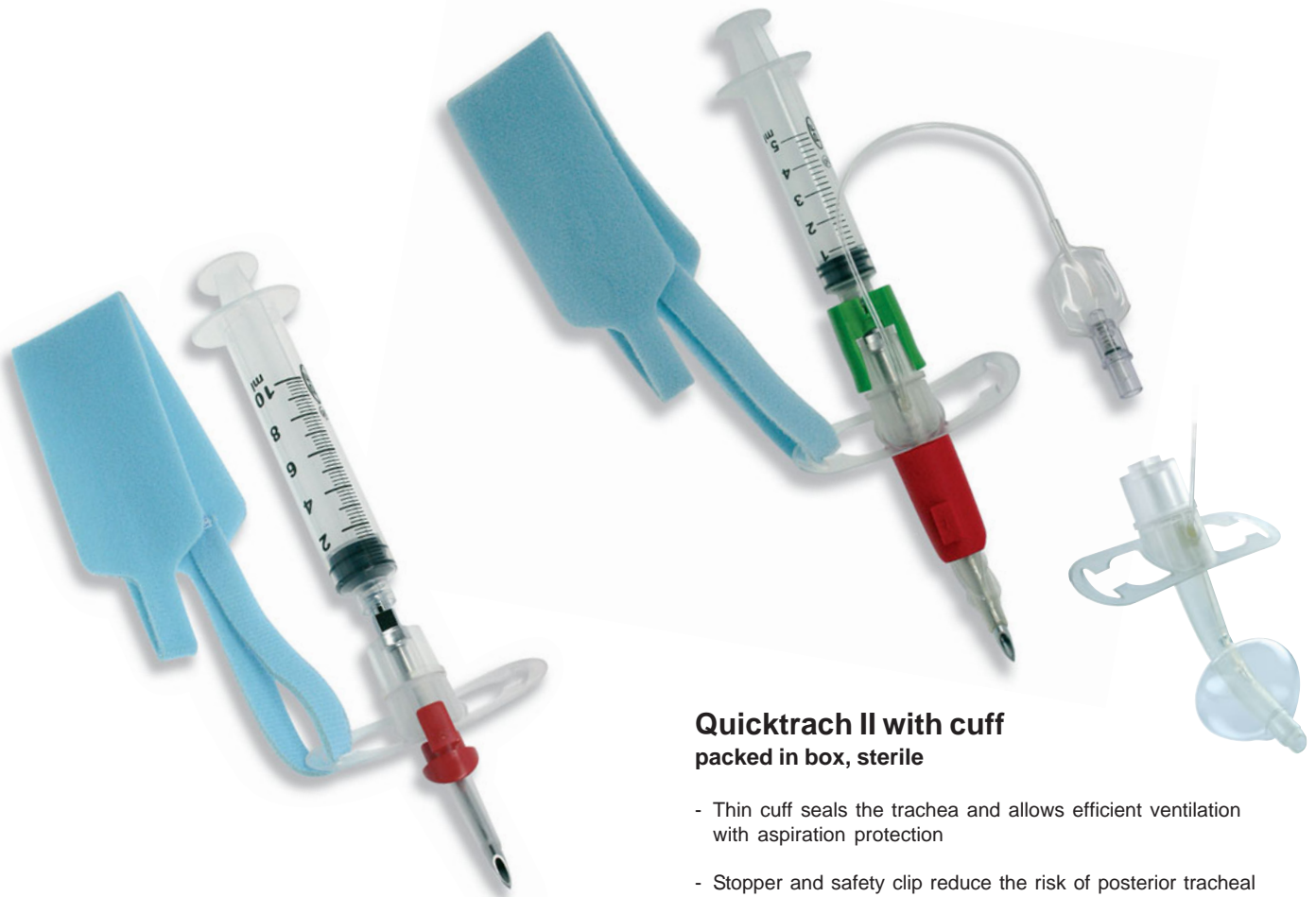

Quicktrach I packed in box, sterile

Patient	I.D.	REF	Box
Infant	1.5 mm	30-04-015-1	1
Child	2.0 mm	30-04-002-1	1
Adult	4.0 mm	30-04-004-1	1

Quicktrach II with cuff packed in box, sterile

- Thin cuff seals the trachea and allows efficient ventilation with aspiration protection
- Stopper and safety clip reduce the risk of posterior tracheal wall perforation
- Anatomically shaped cannula adjusts to the trachea due to „memory effect“

Patient	I.D.	REF	Box
Adult	4.0 mm	30-10-004-1	1

packed in plastic tube, sterile

- ideal for emergency bags (compact and robust)



Quicktrach I

Patient	I.D.	REF	Box
Child	2.0 mm	30-04-902-1	1
Adult	4.0 mm	30-04-904-1	1

Quicktrach II with Cuff

Patient	I.D.	REF	Box
Adult	4.0 mm	30-10-904-1	1